


When design and function come together...you get **FUSION!**

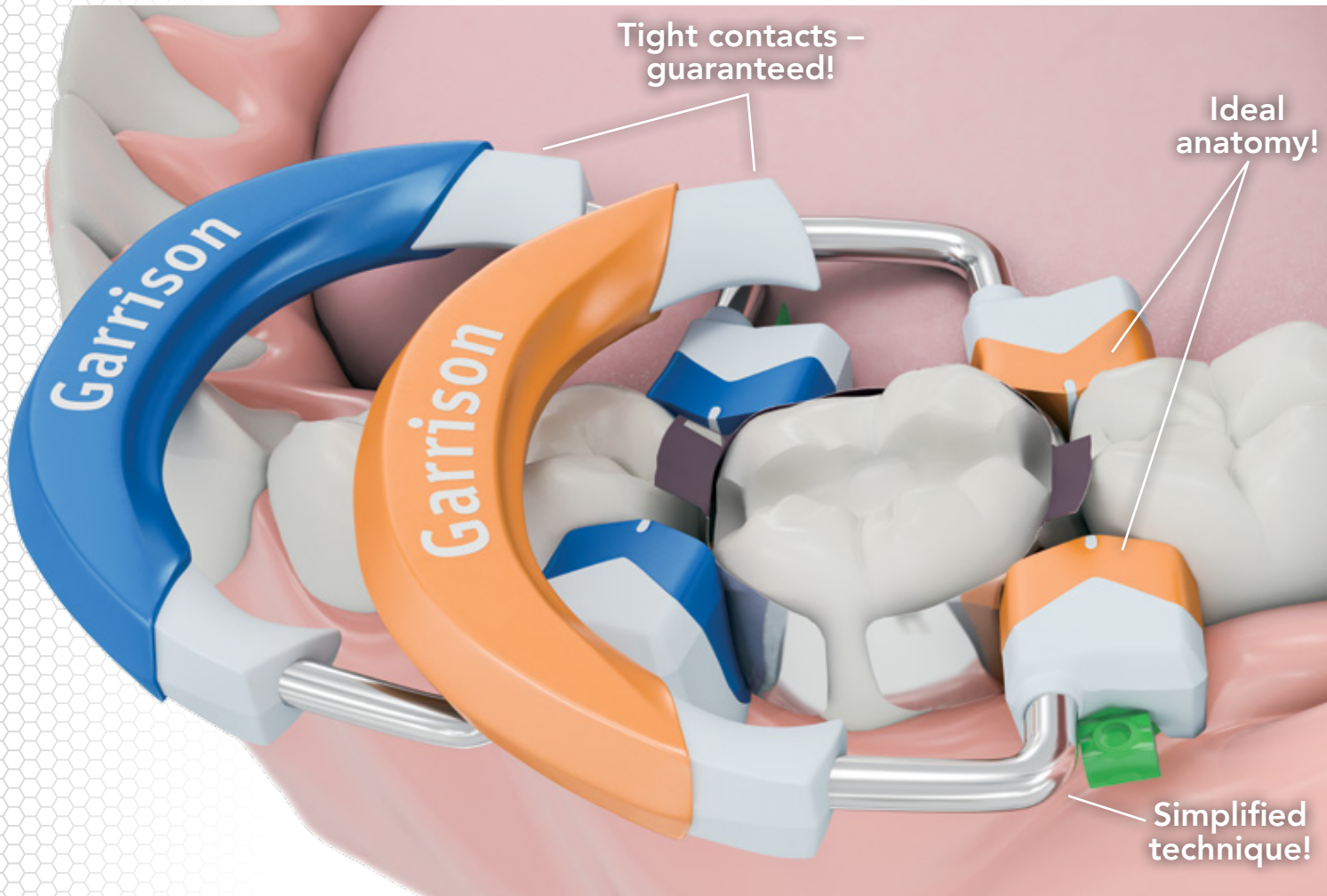

Tight contacts – guaranteed!

Superior tooth separation force through advanced engineering and material selection guarantees perfect contacts every time.

Ideal anatomy!

Soft-Face™ silicone ring tips with enhanced marginal ridge anatomy hug the pre-contoured dead-soft matrices to the tooth for flash free, perfectly contoured restorations.

Simplified technique!


Composi-Tight® 3D Fusion™ is packed with features like easy ring stacking, over-the-wedge placement and Ultra-Grip™ anti-ring slip extensions that make system use intuitive and frustration free.

Composi-Tight® 3D Fusion™ 
Sectional Matrix System

Only a sectional matrix can deliver results like this!

Tofflemire-style System

- ❌ Fails to restore proximal anatomy
- ❌ Thin contact at the marginal ridge
- ❌ Large food trap below
- ❌ Increased likelihood of fracture, occlusal interference, recurrent caries and periodontal disease

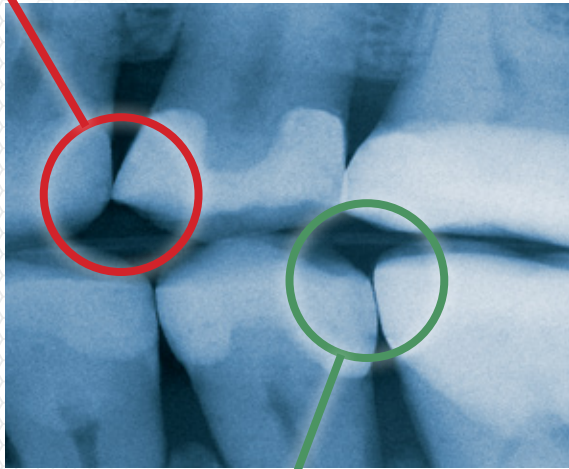

Composi-Tight® Sectional Matrix System

- ✓ Operator-friendly retaining system
- ✓ Naturally contoured bands
- ✓ Anatomically correct contacts
- ✓ Contacts at the height of contour
- ✓ Rings produce optimal tooth separation for excellent, tight contacts



Only Composi-Tight® 3D Fusion™ can make the technique this intuitive and frustration free!

1 Place the matrix band:




Dead-soft three dimensionally contoured matrices have placement tabs for easier occlusal placement and Garrison's revolutionary non-stick coating for hassle-free removal.

2 Wedge firmly:



Composi-Tight® 3D Fusion™ Ultra-Adaptive wedges glide in but won't back out. Their soft, flexible exterior readily adapts to interproximal irregularities for a superior gingival seal.

3 Place the ring:

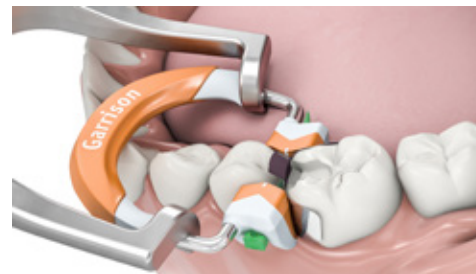

All three of Composi-Tight® 3D Fusion's separator rings will fit right over the wedge for easier placement. With Ultra-Grip™ retention extensions, they will stay firmly in place from the distal of the canine to the widest of molar preparations.

An all-new system designed to handle even the toughest cases

More than 20 years of sectional matrix leadership and an incredible amount of feedback from users all over the globe were poured into the design. All this hard work has resulted in separator rings with significantly improved performance and ease of use.



Get a grip! All Composi-Tight® 3D Fusion™ separator rings incorporate Ultra-Grip™ retention extensions that anchor them securely to virtually any tooth. These rings stay where you place them!



Hugs the curves –
Soft-Face™ silicone tips with enhanced marginal ridge anatomy hug the matrix band to the tooth for flash free restorations




It just keeps going...
Enhanced durability for extended ring life and greater tooth separation



Stack it up – Tru-Stack™ technology allows placement of the orange ring (tall) over the blue ring (short) for easier MOD and multiple tooth restorations


Wider indications, wider smiles!

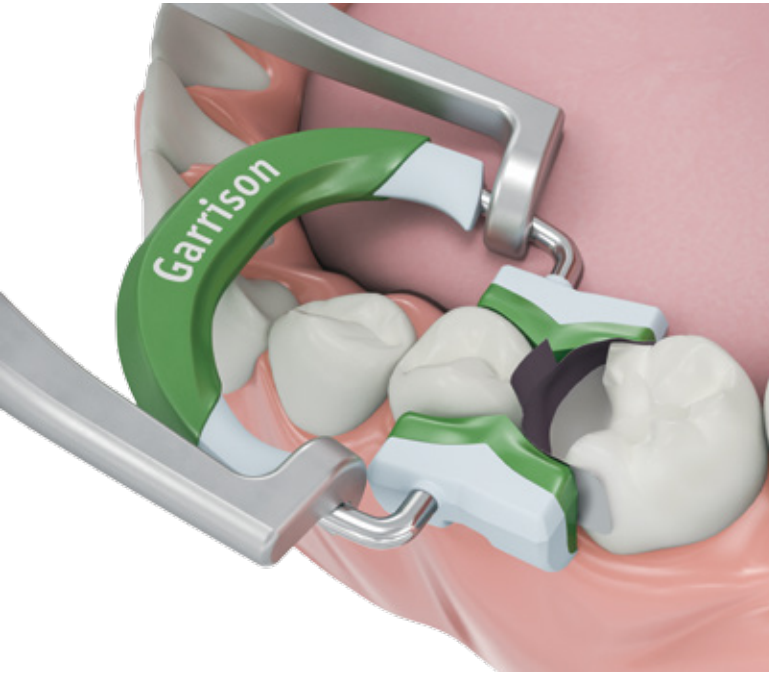


Dentistry's first Wide Preparation separator ring (green) makes Composit-Tight® 3D Fusion™ the most versatile and user-friendly sectional matrix system available. The Wide Prep ring greatly simplifies what until now has been one of the most challenging posterior composite restorations. Now you can restore these big preps with confidence knowing that you will achieve excellent contact and contour without having to take extra steps.

Matrix band perfection

Composit-Tight® 3D Fusion™ Full Curve Non-Stick Matrices are three dimensionally contoured, ultra-thin and dead-soft to help you recreate ideal interproximal anatomy.

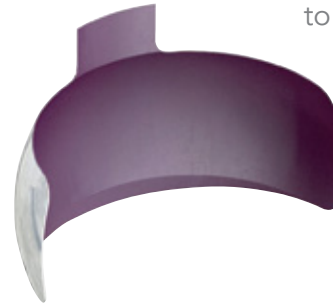

Easy in – easy out! A Grab-Tab™ centered on the occlusal edge of the band simplifies placement while Slick Bands™ non-stick coating reduces composite adhesion by an impressive 92%!

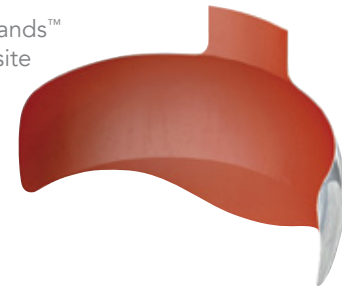


Ultra Adaptive Wedges

Combining Soft-Face™ adaptive materials, a firm inner core and advanced mechanical features to produce a wedge that truly works. Every time.



Fins with serious function
Soft retentive fins smoothly fold down during wedge insertion and then spring back when clear of the interproximal space. Wedge back-out is a thing of the past.

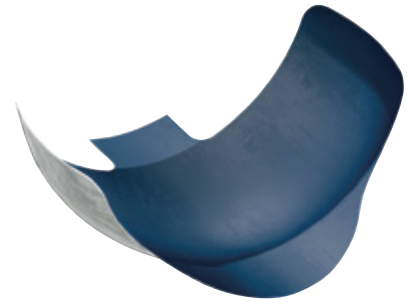

Go deep or wide – Improved subgingival extensions and added wrap-around length make for the perfect matrix band to use with the green Wide Preparation ring. Now you're ready to tackle even the widest preps!



A real softie Compositight® 3D Fusion's Soft-Face™ over-mold allows the wedge to do what no other wedge can truly do – actually adapt to interproximal irregularities.

Tough on the inside
The firm inner skeleton allows for easy insertion and tooth separation just like a traditional wedge.

Compositight® 3D Fusion™ forceps for precise control

Compositight® 3D Fusion™ includes dentistry's premier drop-forged ring placement forceps. Regardless of your hand size or strength, Compositight® 3D Fusion™ forceps will give you incredible control for the precise placement of super-strong separator rings over even the widest molars.



Ergonomic for the masses – Whether your hands are big or small, Compositight® 3D Fusion™ forceps give you the leverage you need for precise ring placement.

Long-term investment – Drop-forged stainless steel prevents leverage-robbing flexion and makes these forceps incredibly durable.

Functionality for all – As an added bonus, Compositight® 3D Fusion™ forceps will work with virtually any ring from any manufacturer.



Only Garrison® can bring you a comprehensive introductory kit like this!

Composi-Tight® 3D Fusion™ Kit FX-KFF-10

Contents:

- 1 – Short Ring (blue)
- 1 – Tall Ring (orange)
- 1 – Wide Prep Ring (green)
- 70 – Assorted Matrix Bands
- 80 – Assorted Ultra-Adaptive Wedges
- 12 – Assorted FenderWedges®
- 1 – Ring Placement Forceps



All Garrison products are covered by our six-month money-back satisfaction guarantee. If you are not satisfied with your product, return the unused portion and we will issue a refund. It's that simple.

Garrison
Dental Solutions

toll free 888.437.0032
www.garrisondental.com